


# MicroMega HERO 642

ROOT CANAL SHAPING

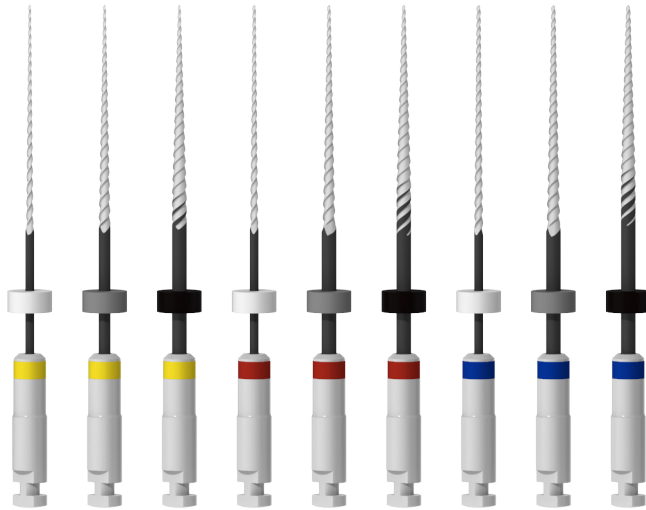

Root canal shaping sequence

# MicroMega

# HERO 642

The HERO 642 is a variable-taper Nickel-Titanium instrument for root canal shaping

- The large «residual core» increases resistance to fracture
- Progressive pitch reduces screwing tendency
- Presence of an inactive tip
- Crown-down technique



Wire diameters

0.60mm

0.80mm

1.00mm



Diameter

20

à

30



Tapers

2%

4%

6%



Available lengths

21mm

25mm



Continuous rotation



Motor speed

300rpm

to

600rpm

# MicroMega Hero 642

6 instruments, non sterile pack

References	020	025	030
P6 HERO 642 L21 2% Classics	20136114	20136115	20136116
P6 HERO 642 L21 4% Classics	20136120	20136121	20136122
P6 HERO 642 L21 6% Classics	20136110	20136111	20136112
P6 HERO 642 L25 2% Classics	20136101	20136102	20136103
P6 HERO 642 L25 4% Classics	20136107	20136108	20136109
P6 HERO 642 L25 6% Classics	20136132	20136133	20136134

## Best-sellers

HERO642 N20 6%



HERO642 N20 4%



HERO642 N20 2%




HERO642 N25 4%



HERO642 N30 4%



# Visit the product page

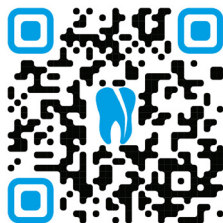

[micro-mega.com](https://micro-mega.com)

# Follow us on socials!



@micromega\_international



Micro-Mega SA



Micro-Mega International