



# MicroMega One RECI Protocol

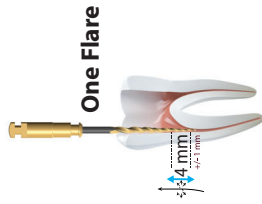


**One Flare**

250-400 rpm  
Max torque 3 N.cm

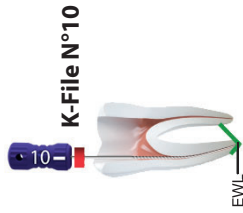

**K-File N°10**

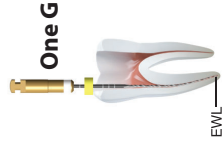


**One G**

250-400 rpm  
Max torque 1.2 N.cm



**K-File N°15**

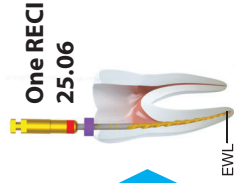

**One RECI  
25.06**

400 rpm  
Max torque 4N.cm  
170° ↺ 60°



Protocol movies  
for all sizes

# MicroMega One RECI Protocol Card

Visit our website for all the latest news:



[micro-mega.com](http://micro-mega.com)

Follow us on socials!



@micromega\_international Micro-Mega SA



Micro-Mega  
International

