


MicroMega One Curve / One Curve mini Protocol



One Flare


250-400 rpm
Max torque 3 N.cm

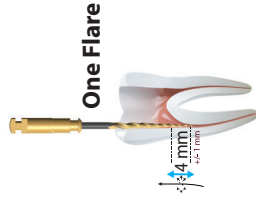

K-File N°10




EWL

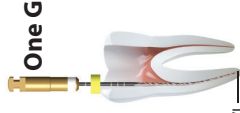

One G

250-400 rpm

Max torque 1.2 N.cm

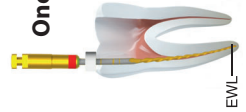


One Curve 25.04

250-400 rpm

Max torque 2.5 N.cm




K-File N°15

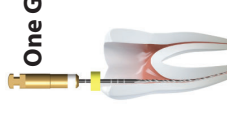

One G

250-400 rpm

Max torque 1.2 N.cm

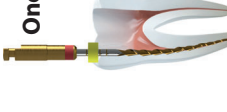

One Curve 25.06

250-400 rpm

Max torque 2.5 N.cm



MicroMega

One Curve / One Curve mini
Protocol Card

Visit our website for all the latest news:



micro-mega.com

Follow us on socials!



@micromega_international | Micro-Mega SA



Micro-Mega
International

